


Checker[®]
handheld colorimeter

ACCURATE
PHOSPHATE
MEASUREMENT

Adaptation of Standard Method
4500-P E

WWW.HANNACHECKER.COM

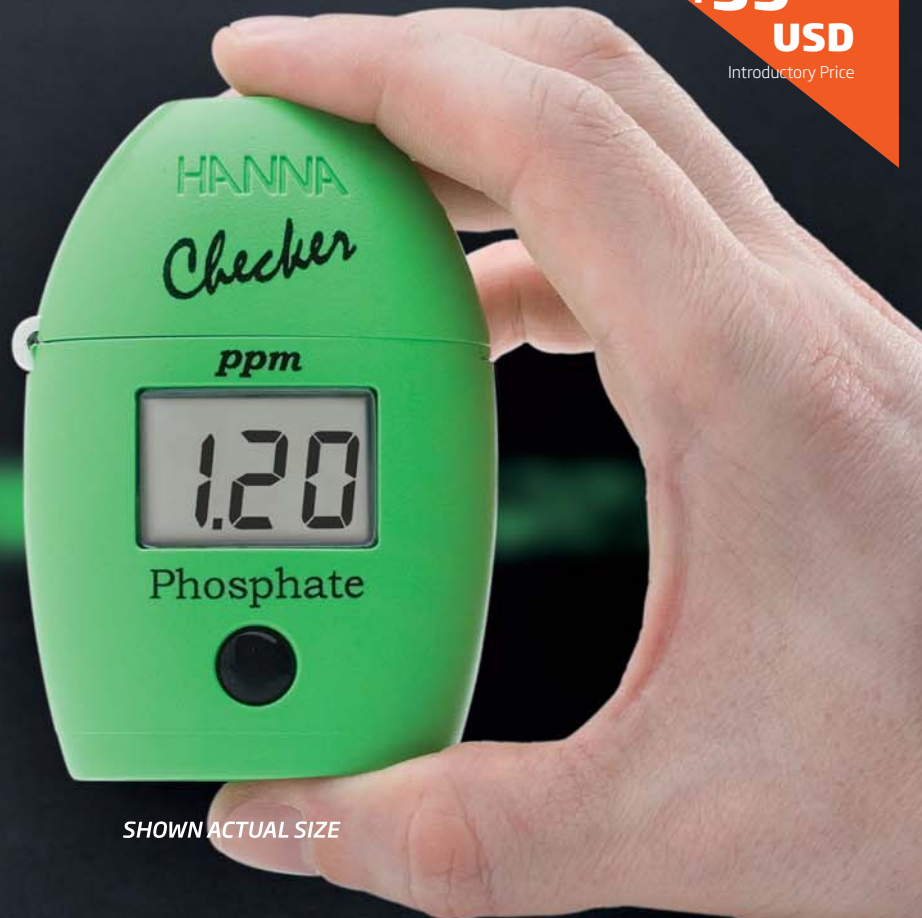

SHOWN ACTUAL SIZE

NEW
\$39.50
USD

Introductory Price

 **HANNA**[®]
instruments
With Great Products, Come Great Results™

Checker[®] HC

handheld colorimeter

Easier to use and more accurate than chemical test kits

- Adaptation of Standard Method 4500-P E
- Accuracy $\pm 4\%$ of reading ± 0.04 ppm (mg/L)
- 0.01 ppm (mg/L) resolution (250 points)
- Large, easy to read digits
- Auto shut-off

Dedicated to a single parameter

- Designed to work with HANNA's powder reagents
- Uses 10 mL glass cuvettes



Small Size, Big Convenience

- Weighing a mere 64 g (2.25 oz.), the Checker[®]HC easily fits into the palm of your hand or pocket
- Use for quick and accurate on the spot analysis
- Single button operation: Zero and Measure
- Operated by a single AAA battery

Ideal for

- Natural, waste and drinking waters



ACCURATE AND AFFORDABLE PHOSPHATE MEASUREMENT

Orthophosphates are found in natural and waste waters. They are commonly added to drinking water as a corrosion inhibitor. The instantaneous analysis of orthophosphate by colorimetric determination provides rapid results using a standard analysis technique.

The HANNA HI 713 Checker[®]HC bridges the gap between simple chemical test kits and professional instrumentation. Chemical test kits are not very accurate and only give 5 to 10 points resolution while professional instrumentation can cost hundreds of dollars and can be time consuming to calibrate and maintain. The HANNA HI 713 Checker[®]HC is accurate and affordable.

The HI 713 Checker[®]HC portable handheld colorimeter features a resolution of 0.01 ppm (250 points) and $\pm 4\%$ of reading ± 0.04 accuracy. The HI 713 Checker[®]HC uses an adaptation of Standard Method 4500-P E.

The contoured style of this Checker[®]HC fits in your palm and pocket perfectly and the large LCD is easy to read. The auto shut-off feature assures the battery life will not be drained if you forget to turn it off.

The HI 713 Checker[®]HC is extremely simple to use. First, zero the instrument with your water sample. Next, add the reagent. Last, Place the vial into the HI 713 Checker[®]HC, press the button and read the results. It's that easy.

WE DESIGN, MANUFACTURE, SUPPLY AND SUPPORT ALL OF OUR PRODUCTS.


THAT'S 360° VALUE

Over the past 30 years, **HANNA** has never failed to design fresh and innovative products. The new Checker®HC's are no exception. We are excited to introduce such high accuracy and ease of use in a small, modern design that fits in the palm of your hand.

With an introductory price of \$39.50, our Checker®HC's are a great value.

HANNA is the largest family-owned manufacturer of analytical instrumentation in the world. The new series of Checker®HC's are manufactured in our European state-of-the-art ISO 9001:2000 production facility. Each Checker®HC is CE compliant to EN 61326-1 and EN 61010-1 standards.

When you buy a **HANNA** product, you're not only buying the best value for your money, but you're also adding the benefit of **HANNA**'s unsurpassed customer service and post-sale technical support.

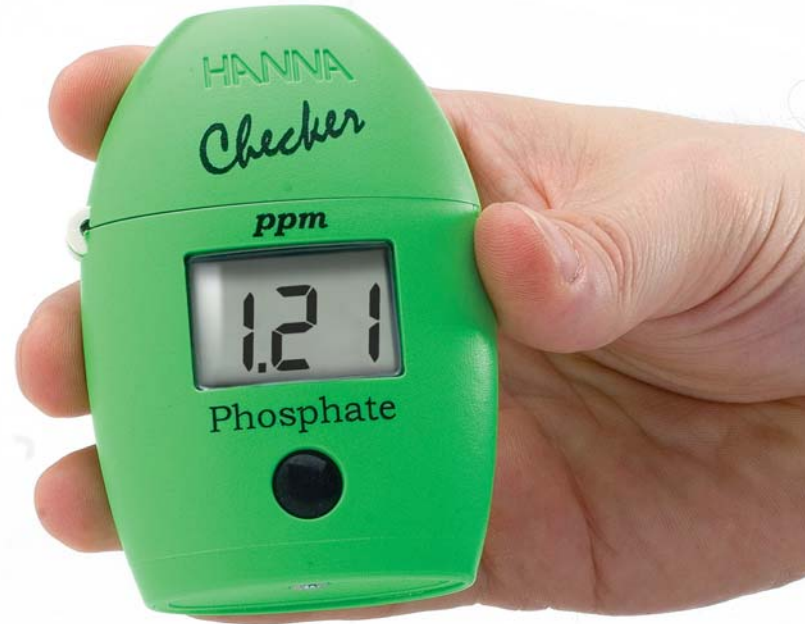

IT'S EASY TO MEASURE SAMPLES WITH THE CHECKER®HC



"Zero" the Checker®HC with your unreacted water sample



Add reagent to your water sample



Place the vial into your Checker®HC



Press the button and read the results. it's that easy!

Checker®
handheld colorimeter

SPECIFICATIONS	HI 713 (Phosphate)
Range	0.00 to 2.50 ppm (mg/L)
Resolution	0.01 ppm (mg/L)
Accuracy	±4% of reading ±0.04 ppm (mg/L) @ 25°C
Light Source	LED @ 525 nm
Light Detector	Silicon photocell
Environment	0 to 50°C (32 to 122°F); RH max 95% non-condensing
Battery Type	(1) 1.5V AAA
Auto-off	After 2 minutes of non-use and after read
Dimensions	81.5 x 61 x 37.5 mm (3.2 x 2.4 x 1.5")
Weight	64 g (2.25 oz.)
Method	Adaptation of Standard Method 4500-P.E.

ORDERING INFORMATION

HI 713 Checker®HC is supplied with (2) sample cuvettes with caps, (6) powder reagents for phosphate, (1) 1.5V AAA battery and instructions.

Total Chlorine (HI 711), **Free Chlorine** (HI 701), **Iodine** (HI 718) and **Iron** (HI 721) Checker®HC's also available

REAGENTS AND STANDARDS

HI 713-25 Reagents for 25 tests (phosphate)


ACCESSORIES

HI 731318 (4) Cuvette cleaning cloth

HI 731331 (4) Glass cuvettes

HI 731353 (4) Caps for cuvettes

HI 93703-50 Cuvette cleaning solution, 230 mL

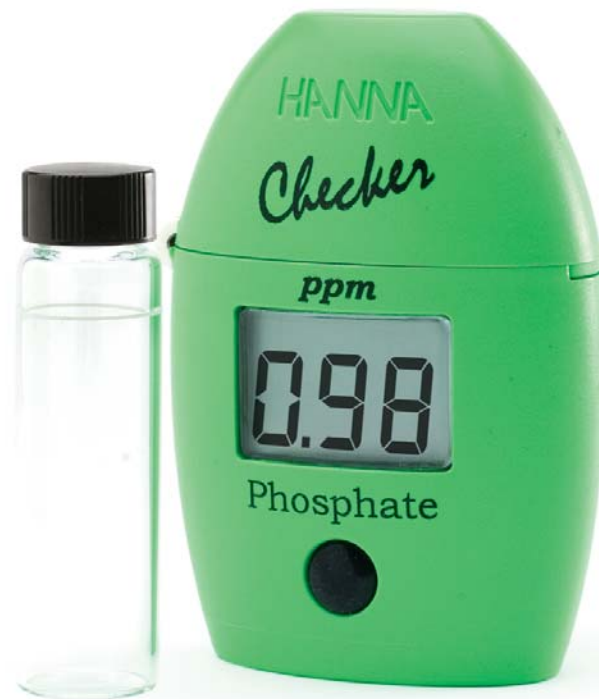

SHOWN ACTUAL SIZE

WWW.HANNACHECKER.COM

HANNA®
instruments
With Great Products, Come Great Results™